

100bp DNA Ladder(Catalog# TBS4030)

DESCRIPTION

100 bp DNA ladder is composed of nine bands of 100 bp, 200 bp, 300 bp, 400 bp, 500 bp, 600 bp, 800 bp, 1,000 bp, and 1,500 bp. The 500-bp band with higher concentration is easily distinguishable from the others.

APPLICATIONS

- DNA marker

KIT CONTENTS


TBS4030-0.5: 0.5 mL x 1
TBS4030-1.0: 0.5 mL x 2

0.5 mL in 10 mM Tris-HCl (pH=7.5), 1 mM EDTA.
Concentration: 250 ng/μL

STORAGE CONDITIONS

For long term storage, keep kit at -20°C.

Example

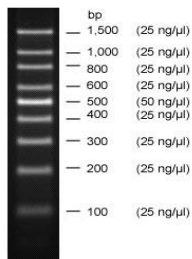

2% agarose gel

RECOMMENDED LOADING VOLUME

5 μL (or 1μL/mm width of gel well) is recommended for use.

NOTE

- 6x DNA loading buffer can be ordered separately.

QUALITY CONTROL

- Agarose gel analysis shows that the bands between 100 bp and 1,500 bp are accurate in size and distinguishable. The band of 500 bp is more easily identified than any other bands.
- 5 μL quantities of the DNA marker do not show any visible degradation after incubation in restriction enzyme buffer overnight at 37°C.

RELATED PRODUCTS

6x or 10x DNA Loading Buffer (Catalog#TBS5018 or TBS5023)
Tail DNA Extraction kit (Catalog# TBS6005)
Fast DNA Extraction kit (Catalog# TBS6008)
2x Regular PCR Kit (Catalog#TBS4004)
2x PCR Red Mix kit (Catalog# TBS4005)

This product is for *in vitro* research use only and is not intended for use in humans or animals in therapeutic or diagnostic procedures.